## BACINELLA SHAMPOO SHAMPOO BASIN

### Manuale d'uso - User manual

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**ATTENTION:** The operators must carefully read and completely understand the present manual before using the product.

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## INTENDED USE

1. Designed for bedridden or immobilized persons.

## CONTENTS

1. Shampoo basin with drainage hose attached.
2. Water pouch with shower head and clip.
3. A plastic pipe and string for better water pouch handling and hanging.
4. A repair kit.

1. Shampoo basin

2. Water pouch kit

3. A plastic pipe and string

3. Place the drainage hose to a water bucket (not included).


**NOTE:** water bucket must be placed lower than the shampoo basin to ensure soiled water flows downwards.

4. Unplug the drainage plug for water drainage while in use.

Drainage plug

5. Place shampoo basin under patient's neck area, adjust for comfort when necessary.

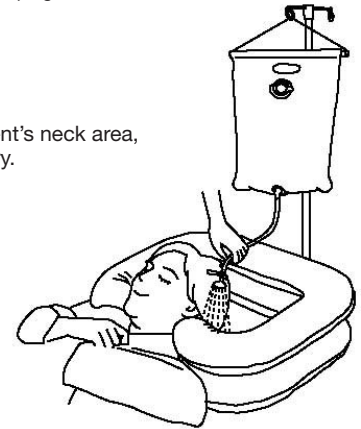Application illustration

## FEATURES

1. A 20cm deep shampoo basin with dual tubing for sturdiness and comfort.
2. It comes with a drainage plug with hose for ease of discharging soiled water during or after use.
3. A separate water pouch for easy water transport where water source may not be conveniently accessible. It also comes with a shower head and an adjustable clip for ease of water flow control.

## SETUP INSTRUCTION

1. Inflate the shampoo basin either by mouth or with air pump (not included). **NOTE:** the air valve has a one-way valve inside to help to prevent air from leaking out. To inflate, squeeze air valve by its neck to help to open the one-way valve inside.
2. Place the shampoo basin on a flat and semi-rigid surface for better use.

Airflo (Xiamen) Medical Co., Ltd.  
4F, No.6, East Haijian Road, Haicang, Xiamen, Fujian, China  
Made in China

Sungo Certification Co., Ltd.  
RM101, MAPLE HOUSE, 118 HIGH STREET, PURLEY, LONDON, ENGLAND