
GIMA

PROFESSIONAL MEDICAL PRODUCTS

Gima S.p.A. - Via Marconi, 1 - 20060 Gessate (MI) Italy
Italia: tel. 199 400 401 - fax 199 400 403
Export: tel. +39 02 953854209/221/225 fax +39 02 95380056
gima@gimaitaly.com - export@gimaitaly.com
www.gimaitaly.com

BILANCIA GIMA PROFESSIONAL DIGITALE
GIMA PROFESSIONAL DIGITAL SCALE
BALANCE NUMÉRIQUE PROFESSIONAL GIMA
DIGITALWAAGE „PROFESSIONAL“ GIMA
BÁSCULA PROFESIONAL DIGITAL GIMA
BALANÇA PROFISSIONAL DIGITAL GIMA
ΨΗΦΙΑΚΗ ΕΠΑΓΓΕΛΜΑΤΙΚΗ ΖΥΓΑΡΙΑ GIMA
(GIMA) "ميزان مهني رقمي " جيما "

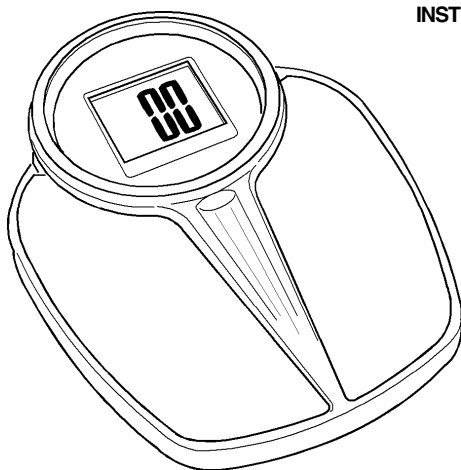MANUALE D'USO E MANUTENZIONE
USE AND MAINTENANCE BOOK
INSTRUCTIONS DE FONCTIONNEMENT ET ENTRETIEN
BETRIEBS UND WARTUNGS ANWEISUNGEN
MANUAL DE USO Y MANTENIMIENTO
MANUAL DE USO E MANUTENÇÃO
ΕΓΧΕΙΡΙΔΙΟ ΧΡΗΣΗΣ ΚΑΙ ΣΥΝΤΗΡΗΣΗΣ
دليل الإستعمال والرعاية

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FEATURES

The GIMA Professional digital scale is a device for measuring the weight of the human body. The wide measurement range makes it ideal for any kind of use.

Technical Features

Display LCD 2,5"

Maximum capacity: 200 kg / 440 lb

Precision: 100 g / 0,2 lb

Permitted tolerance: <100 kg ±0,5 kg >100 kg ±0,7 kg

Automatic switch off

Operated with 2 batteries type CR 2032

Not certified product class III

PRESCRIPTIONS

Do not use the equipment in case it is damaged. Apply to your retailer.

Avoid precarious repairs. Repairs shall be carried out with original spare parts only, which shall be installed according to the intended use.

Since the product is made of corrosion-proof materials suitable for the environmental conditions foreseen for its normal use, does not require special care, however it is necessary to store it in a closed place making sure that is protected from dust and dirt to assure its hygienic conditions. Moreover, it is recommended to store the product in a place which can be reached easily by the personnel in case of necessity.

UNPACKING

Always remember that packing elements (paper, cellophane, stitches, adhesive tape, etc.) can cut and/or hurt if they are not carefully handled. They shall be removed with adequate means and shall not be left at the mercy of irresponsible persons; the same is valid for tools used to remove packages (scissors, knives, etc.).

After opening the packages, first of all it is necessary to check all pieces and parts composing the product. Check that they are all present and in perfect conditions.

Before using the scale remove the isolation tab from batteries.

OPERATION

- Select the unit of measure (kg, lb or st) using the selector on the bottom of the scale
- Position the scale on the floor in a stable position (ex.: not on a rug)
- Switch on the scale by touching it with your foot
- Wait approx. 3 seconds until the display shows "0.0". At this point the scale is ready for use.
- Step onto the scale and verify that both feet are positioned correctly
- Try to stand as still as possible for a precise measurement
- The measured value will flash on the display twice and will remain on the display for approx. 10 seconds
- After stepping down from the scale it will automatically switch off after approx. 6 seconds
- To obtain valid results for statistical purposes, weigh yourself every day possibly at the same time. It is best to weigh yourself without clothes and in the morning as soon as you wake up.

Caution: If the weight on the scale exceeds 200 kg the display will show "ERR". Do not overload the scale or it could get damaged!

Do not use the scale near emissions of electromagnetic waves (mobile phones, microwave ovens).

MAINTENANCE

The device was designed to last a lifetime. It is resistant to chemical attack, except for M.E.K. ketons, diluting solutions and polish removers.

When the display shows "LO" replace the batteries (2 x CR 2032). The battery box is on the bottom of the scale.

Disposal: *The product must not be disposed of along with other domestic waste. The users must dispose of this equipment by bringing it to a specific recycling point for electric and electronic equipment.*

For further information on recycling points contact the local authorities, the local recycling center or the shop where the product was purchased. If the equipment is not disposed of correctly, fines or penalties may be applied in accordance with the national legislation and regulations.

GIMA WARRANTY CONDITIONS

Congratulations for purchasing a GIMA product.

This product meets high qualitative standards both as regards the material and the production.

The warranty is valid for 12 months from the date of supply of GIMA.

During the period of validity of the warranty, GIMA will repair and/or replace free of charge all the defected parts due to production reasons. Labor costs and personnel traveling expenses and packaging not included.

All components subject to wear are not included in the warranty.

The repair or replacement performed during the warranty period shall not extend the warranty.

The warranty is void in the following cases: repairs performed by unauthorized personnel or with non-original spare parts, defects caused by negligence or incorrect use.

GIMA cannot be held responsible for malfunctioning on electronic devices or software due to outside agents such as: voltage changes, electro-magnetic fields, radio interferences, etc.

The warranty is void if the above regulations are not observed and if the serial code (if available) has been removed, cancelled or changed.

The defected products must be returned only to the dealer the product was purchased from.

Products sent to GIMA will be rejected.

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