BILANCIA BODY FAT LIBRA
LIBRA BODY FAT SCALE
PÈSE-PERSONNES BODY FAT LIBRA
BÁSCULA BODY FAT LIBRA
BALANÇA BODY FAT LIBRA
LIBRA WAGA BODY FAT

MANUALE D'USO E MANUTENZIONE
USE AND MAINTENANCE BOOK
INSTRUCTIONS DE FONCTIONNEMENT ET ENTRETIEN
MANUAL DE USO Y MANTENIMIENTO
MANUAL DE USO E MANUTENÇÃO
PODRECZNIK EKSPLOATACJI I KONSERWACJI

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Features

- * Capacity: 150 kg / 330 lb
- * Graduation: 0.1 kg / 0.2 lb
- * Auto on or press Unit to switch on
- * 3 function key: SET, ▲, ▼.
- * It may measure body fat, water, muscle and bone
- * Identify automatically function: can store 8 users' parameter and the last measure data
- * Auto on and Auto off function
- * Over load indication "-----"
- * Power: 2x CR 2032 Batteries

Display LCD

To change batteries

- Ensure +/- terminals are correct.
- Use proper tool to help release the battery from the compartment if necessary.

To input personal information

- * Set up personal data(gender/age/stature) before using, it can set up 8 users' data

1. Press "set" to enter the mode of setup, the flashing Pn (n can be 1~8) means the memory of the data. Press ▲ & ▼ to change number, gender, stature, age, then press "SET" to confirm. The scale will back to the mode of weighing after finishing setting up.

2. Set up the gender: press UP, DOWN (can select female, male, athlete mode), then press SET to confirm

3. Set up the stature: Press UP, DOWN, and then press SET to confirm


4. Set up age: Press UP, DOWN, and then press SET to confirm

To weigh and measure

1. After setting up of the data, when LCD displays 0.0, you may weigh or stand on the scale platform directly, the weight value will be shown after a few seconds.

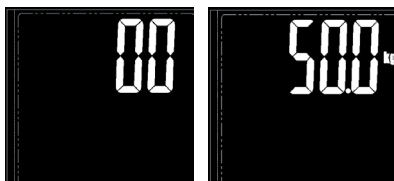2. After the weight is stable, the flashing "0" means that the scale is measuring the fat, the "0" signal moves from left to right for two times, and then shows the data

Display the measurement


After the measuring, if the user don't choose the user No, the scale will automatically identify the users according to the weight. And the LCD will display the weight, fat, water, muscle and bone data for three times.

Automatically identify the user measurement

When the user stand on the scale directly, the scale will turn on and lock the weight.

When the weight is stable, display moving "0" means it is measuring. If it can identify the user successfully, it will display the measurement. If it detects many similar weight, it will display the relative user No., Use up to choose the left user No. and use down to choose the right user No. After finished choosing the user, it will display the measurement and turn off.

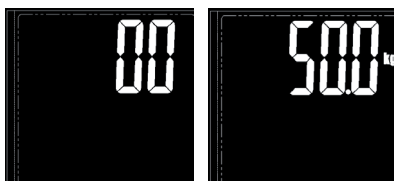If you don't choose the user, it will display the user No. several seconds and turn off.

If it occurred error or it can't identify the user, it will display the present weight and turn off.

Body fat, water and muscle rate range

The body fat ratio (%): 3%-50%

Age	Female				Male			
	Underfat	Healthy	Overfat	Obese	Underfat	Healthy	Overfat	Obese
<39	<21	21.1-33	33.1-39	>39	<8	8.1-20	20.1-25	>25
40-59	<23	23.1-34	34.1-40	>40	<11	11.1-22	22.1-28	>28
>59	<24	24.1-36	36.1-42	>42	<13	13.1-25	25.1-30	>30

The body water ratio(%): 25%-75%

Gender	Female			Male		
Age	Low	Normal	High	Low	Normal	High
10-15	<57	57.1-67	>67	<58	58.1-72	>72
16-30	<47	47.1-57	>57	<53	53.1-67	>67
31-60	<42	42.1-52	>52	<47	47.1-61	>61
61-80	<37	37.1-47	>47	<42	42.1-56	>56

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Disposal: The product must not be disposed of along with other domestic waste. The users must dispose of this equipment by bringing it to a specific recycling point for electric and electronic equipment. For further information on recycling points contact the local authorities, the local recycling center or the shop where the product was purchased. If the equipment is not disposed of correctly, fines or penalties may be applied in accordance with the national legislation and regulations.

GIMA WARRANTY CONDITIONS

Congratulations for purchasing a GIMA product.

This product meets high qualitative standards both as regards the material and the production. The warranty is valid for 12 months from the date of supply of GIMA.

During the period of validity of the warranty, GIMA will repair and/or replace free of charge all the defected parts due to production reasons. Labor costs and personnel traveling expenses and packaging not included.

All components subject to wear are not included in the warranty.

The repair or replacement performed during the warranty period shall not extend the warranty.

The warranty is void in the following cases: repairs performed by unauthorized personnel or with non-original spare parts, defects caused by negligence or incorrect use.

GIMA cannot be held responsible for malfunctioning on electronic devices or software due to outside agents such as: voltage changes, electro-magnetic fields, radio interferences, etc.

The warranty is void if the above regulations are not observed and if the serial code (if available) has been removed, cancelled or changed.

The defected products must be returned only to the dealer the product was purchased from. Products sent to GIMA will be rejected.

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