




GVS Filter Technology UK Ltd.
NFC House, Vickers Industrial Estate
Mellishaw Lane
Morecambe, Lancashire
LA3 3EN
United Kingdom
Tel: +44 (0) 1524 847600
Email: gvsuk@gvs.com

FILTER TECHNOLOGY



Instructions For Use; Mode d'emploi; Gebrauchsanweisung; Istruzioni per l'uso; Instrucciones de uso; Gebruiksaanwijzing; Instruções de Utilização; Bruksanvisning; Bruksanvisning; Bruksanvisning; Käyttöohjeet; Instrukcja użytkowania; Инструкция за употреба; Upute za upotrebu; Návod k použití; Kasutusjuhend; Οδηγίες χρήσης; Használati útmutató; Treoracha Úsáide; Lietošanas instrukcija; Naudojimo instrukcijos; Istruzzjonijiet dwar l-Użu; Instrucțiuni de utilizare; Návod na použitie; Navodila za uporabo.

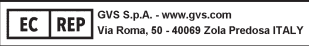

CE 2797

	UK	FRANCE	GERMANY	SPAIN	NETHERLANDS	USA
Classification	Class IIa according to EC Directive 93/42/EEC Annex IX	Classe IIa selon la Directive CE 93/42/CEE Annexe IX	Klasse IIa gemäß EG-Richtlinie 93/42/EWG Anhang IX	Clase IIa según la Directiva 93/42/CEE, Anexo IX	Klasse IIa volgens EU-richtlijn 93/42/EEG Appendix IX	Class II in accordance to QSR 820.00
UMDNS code	11-710 Universal Medical Device Nomenclature System	11-710 Système universel de nomenclature des dispositifs médicaux	11-710 Nomenklatur für Medizingeräte	11-710 Sistema universal de nomenclatura de dispositivos médicos	11-710 Nomenclatuur voor medische apparaten	11-710 Universal Medical Device Nomenclature System
Packaging code	"XXXX" defines packaging variation e.g. ABUA	"XXXX" définit les emballages variation par exemple ABUA	"XXXX" definiert Verpackung Variation z.B. ABUA	"XXXX" define embalaje variación por ejemplo, ABUA	"XXXX" definitieert verpakking variatie bijv. ABUA	"XXXX" defines packaging variation i.e. ABUA

REF Product Code; Code du produit; Artikelnummer; Codice prodotto; Código de producto; Productcode; Código de produto; Produktkod; Produktkode; Produktkode; Tuotekoodi; Kod produktu; Код на продукта; Sifra proizvoda; Kód výrobku; Toote kood; Kωδικός προϊόντος; Termékkód; Cód Táiige; Izstrādājuma kods; Gaminio kodas; Kódići tal-prodott; Cod produs; Код výroку; Koda izdelka.	2000/53XXXX Suction Filter	2200/02XXXX Suction Filter	2200/02DIKXXXX Suction Set	2200/06XXXX Suction Filter	2200/16XXXX Suction Filter	2200/21XXXX Suction Filter	2200/26XXXX Suction Filter	2200/35XXXX Suction Filter	2200/36XXXX Suction Filter 3A1385	2200/55XXXX Suction Filter	2200/70XXXX Suction Filter	2200/902XXXX 90mm Suction Filter	2200/911XXXX 90mm Suction Filter 3A1925
Product description; Description du produit; Produktbeschreibung; Descripción producto; Descripción de producto; Productomschrijving; Descrição do produto; Produktbeskrivning; Produktbeskrivelse; Produktbeskrivelse; Tuotteen kuvaus; Opis produktu; Описове на продукта; Opis proizvoda; Opis výrobku; Toote kirjeldus; Тарыповор; Тарыповор; Termékleírás; Tuarisc ar Tháirge; Izstrādājuma apraksts; Gaminio aprašymas; Deskripcijoni tal-prodott; Descriera produsului; Opis výrobku; Opis izdelka. Suction Filters and Sets; Filtrés et ensembles d'aspiration; Saugfilter and Sets; Filtri e set di aspirazione; Filtros y conjuntos de succión; Aanzuigfilters en sets; Filtros e conjuntos de sucção; Sugfilter och uppsättningar; Sugefiltre og -sæt; Sugefiltre og sett; Imusuodattimet ja -sarjat; Filtry i zestawy ssące; Всмокателни филтри и комплекти; Usisni filtri i pribori; Sacci filtry a sady; Imemisfiltrid ja komplekti; Φίλτρα και σετ αναρρόφησης; Szívószűrők és -készletek; Foirne agus Scagairi Súcáin; Iesúksanas filtri un komplekti; Siurbimo filtral ir rinkiniai; Filtri u Setti-jiet tal-Gbid tal-Arja; Filtre și seturi de aspirare; Sacie filtre a súpravy; Sesalni filtri in kompleti.	Inlet/Outlet Connectors: 8mm-11mm hose barbed both sides, machine side indicated by markings on outer ring. Approximate dimensions: 54mm diameter x 51mm height. Flow direction indicated by an arrow on the yellow overmould ring of the filter. Weight: 11g (approx.). Unidirectional filter.	Inlet/Outlet Connectors: OD 8mm (approx.) hose barbed and ID 4mm (approx.) - both sides. Approx. dimensions 65mm diameter x 51.5mm length. Flow direction indicated by IN marking on filter face. Weight 17gm (approx.). Unidirectional filter.	Inlet/Outlet Connectors: OD 8mm (approx.) hose barbed and ID 4mm (approx.) both sides with luer lock connector. Approx. dimensions of 2200/01: 65mm diameter x 51.5mm length. Approx. dimensions of tube: 30mm length, 9mm OD x 6mm ID. Flow direction indicated by IN marking on filter face. Weight of filter: 17gm (approx.).	Inlet/Outlet Connectors: 5-9.5mm hose barbed - both sides. Approx. dimensions 65mm diameter x 65.8mm height. Flow direction indicated by IN marking on filter face. Weight 17gm (approx.). Unidirectional filter.	Inlet/Outlet Connectors: 11mm to 15mm hose barbed both sides. Approx. dimensions 64.7mm diameter x 55.7mm height. Flow direction indicated by IN marking on filter face. Weight 17gm (approx.). Unidirectional filter.	Inlet/Outlet Connectors: IN 11mm -15mm hose barbed, OUT 15mm Male ISO. Approx. dimensions 65mm diameter x 56mm length. Flow direction indicated by IN marking on filter face. Weight 17gm (approx.). Unidirectional filter.	Inlet/Outlet Connectors: IN 8mm hose barbed, OUT 15mm Male ISO. Approx. dimensions 65mm diameter x 56mm length. Flow direction indicated by IN marking on filter face. Weight 17gm (approx.). Unidirectional filter.	Inlet/Outlet Connectors: IN 8mm hose barbed, OUT 7mm extended hose taper. Approx. dimensions 65mm diameter x 63.4mm length. Flow direction indicated by IN marking on filter face. Weight 17gm (approx.). Unidirectional filter.	Inlet/Outlet Connectors: IN 7mm hose plain, OUT 8mm hose barbed. Approx. dimensions 65mm diameter x 52mm length. Flow direction indicated by IN marking on filter face. Weight 17gm (approx.). Unidirectional filter.	Inlet/Outlet Connectors: IIN 8.6mm hose plain, OUT 12mm hose straight. Approx. dimensions 65mm diameter x 42.2mm height. Flow direction indicated by IN marking on filter face. Weight 17gm (approx.). Unidirectional filter.	Inlet/Outlet Connectors: IN 8.1mm conical, OUT 12mm hose barbed. Approx. dimensions 65mm diameter x 60.3mm height. Flow direction indicated by IN marking on filter face. Weight 17gm (approx.). Unidirectional filter.	Inlet/Outlet Connectors: 11mm to 15mm hose barbed both sides. Approx. dimensions 92mm diameter x 61mm height. Flow direction indicated by IN marking on filter face. Weight 47gm (approx.). Unidirectional filter.	Inlet/Outlet Connectors: IN 11mm hose plain, OUT 11-15mm hose barbed. Approx. dimensions 92mm diameter x 61.8mm height. Flow direction indicated by IN marking on filter face. Weight 47gm (approx.). Unidirectional filter.
Product supplied non sterile; Produit livré non stérile; Das Produkt wird nicht steril geliefert; Prodotto non sterile; Producto suministrado sin esterilizar; Product niet-steriel aangeleverd; Produto fornecido não-estéril; Produkten levereras ej steril; Produkt leveret ikke-steril; Produkt leverit ikke-steril; Tuote toimitetaan ei-steriilinä; Produkt dostarczany w stanie niesterilnym; Προϊόνταται σε άστυαν ηεστερηνη; Προϊονδ isporuđen nesterilnim; Výrobek je dodávan v nesterilnim stavu; Toode tarnitakse mittesteriilsena; Το προϊόν παράγεται μη στυείο; A termék nem sterilen szállitva; Soláthrailear an táirge neamhshteiriúil; Izstrādājums piegādāts nesterilis; Gaminys tiekiamas nesterilus; Produkt format mhux steril; Produs furnizat nesteril; Výrobok sa dodáva v nesterilnom stave; Izdelek je dobavljen nesteriliziran.	Yes; Oui; Ja; Si; Si; Ja; Sim; Ja; Ja; Ja; Joo; Tak; DA; Da; Ano; Jah; Nai; Igen; Tá; Já; Taip; Iva; Da; Ano; Da.	Yes; Oui; Ja; Si; Si; Ja; Sim; Ja; Ja; Ja; Joo; Tak; DA; Da; Ano; Jah; Nai; Igen; Tá; Já; Taip; Iva; Da; Ano; Da.	Yes; Oui; Ja; Si; Si; Ja; Sim; Ja; Ja; Ja; Joo; Tak; DA; Da; Ano; Jah; Nai; Igen; Tá; Já; Taip; Iva; Da; Ano; Da.	Yes; Oui; Ja; Si; Si; Ja; Sim; Ja; Ja; Ja; Joo; Tak; DA; Da; Ano; Jah; Nai; Igen; Tá; Já; Taip; Iva; Da; Ano; Da.	Yes; Oui; Ja; Si; Si; Ja; Sim; Ja; Ja; Ja; Joo; Tak; DA; Da; Ano; Jah; Nai; Igen; Tá; Já; Taip; Iva; Da; Ano; Da.	Yes; Oui; Ja; Si; Si; Ja; Sim; Ja; Ja; Ja; Joo; Tak; DA; Da; Ano; Jah; Nai; Igen; Tá; Já; Taip; Iva; Da; Ano; Da.	Yes; Oui; Ja; Si; Si; Ja; Sim; Ja; Ja; Ja; Joo; Tak; DA; Da; Ano; Jah; Nai; Igen; Tá; Já; Taip; Iva; Da; Ano; Da.	Yes; Oui; Ja; Si; Si; Ja; Sim; Ja; Ja; Ja; Joo; Tak; DA; Da; Ano; Jah; Nai; Igen; Tá; Já; Taip; Iva; Da; Ano; Da.	Yes; Oui; Ja; Si; Si; Ja; Sim; Ja; Ja; Ja; Joo; Tak; DA; Da; Ano; Jah; Nai; Igen; Tá; Já; Taip; Iva; Da; Ano; Da.	Yes; Oui; Ja; Si; Si; Ja; Sim; Ja; Ja; Ja; Joo; Tak; DA; Da; Ano; Jah; Nai; Igen; Tá; Já; Taip; Iva; Da; Ano; Da.	Yes; Oui; Ja; Si; Si; Ja; Sim; Ja; Ja; Ja; Joo; Tak; DA; Da; Ano; Jah; Nai; Igen; Tá; Já; Taip; Iva; Da; Ano; Da.	Yes; Oui; Ja; Si; Si; Ja; Sim; Ja; Ja; Ja; Joo; Tak; DA; Da; Ano; Jah; Nai; Igen; Tá; Já; Taip; Iva; Da; Ano; Da.	Yes; Oui; Ja; Si; Si; Ja; Sim; Ja; Ja; Ja; Joo; Tak; DA; Da; Ano; Jah; Nai; Igen; Tá; Já; Taip; Iva; Da; Ano; Da.
Ambient Conditions: Conditions ambiantes : Umgebungsbedingungen : Condizioni ambientali : Condiciones ambientales : Condições ambientais : Omgivningsförhållanden : Omgivelsesbetingelser : Krav til omgivelsene : Ympäristöolosuhteet : Warunki otoczenia: Условия на околната среда: Uvjeti okoline: Okolní podmínky: Keskkonnatingimused: Συνθήκες περιβάλλοντος; Könyezeti feltételek: Coiinnilacha Comhthimpeallacha: Apkàrtèje apstàkji: Aplinkos sąlygos: Kundizjonijiet Ambjentali: Conditii ambientale: Okolité podmienky: Okolje.													
Temperature during storage; Température de stockage; Temperatur bei Lagerung; Temperatura durante l'immagazzinamento; Temperatura durante el almacenamiento; Temperatur tijdens opslag; Temperatura em armazenamento; Temperatur under förvaring; Temperatur under opbevaring; Temperatur under opbevaring; Lämpötila varastoinnin aikana; Temperatura podczas składowania; Температура по време на съхранение; Temperatura tijekom skladištenja; Teplota během skladování; Temperatur iadustiamise ajal; Ορυκουσία κατά την αποθήκευση; Hőmérséklet tárolás közben; Teocht le linn na stórála; Temperatura uzglabášanas laikà; Temperatura sandėliavimo metu; Temperatura waqt il-hażna; Temperatura în timpul depozitării; Teplota počas skladovania; Temperatura med skladiščenjem.	5 to 40 °C (41-104 °F)	5 to 40 °C (41-104 °F)	5 to 40 °C (41-104 °F)	5 to 40 °C (41-104 °F)	5 to 40 °C (41-104 °F)	5 to 40 °C (41-104 °F)	5 to 40 °C (41-104 °F)	5 to 40 °C (41-104 °F)	5 to 40 °C (41-104 °F)	5 to 40 °C (41-104 °F)	5 to 40 °C (41-104 °F)	5 to 40 °C (41-104 °F)	5 to 40 °C (41-104 °F)
Performance Data: Caractéristiques ventilatoires : Leistungsdaten : Dati di rendimento : Características de rendimiento : Karakteristieke waarden : Dados de desempenho : Prestandauppgifter : Ydelse : Kapasitetsdata : Suorituskykytiedot : Parametry sprężu; Дани на ефективността; Podaci o izvedbi; Údaje o výkonu; Toimivusandmed; Δεδομένα απόδοσης; Teljesítmény adat; Sonrai Feidhmiochta; Veiktspējas dati; Pasirodymo data; Dejlja tal-Prestazzjoni; Date de performanță; Údaje o výkonosti; Podatki o uspešnosti.													
Air Flow; Flux d'air; Luftstrom; Flusso d'aria; Flujo de aire; Luchtstroom; Fluxo de ar; Luftlöde; Luftstrom; Luftstrom; Ilmavirta; Przepływ powietrza; Вдздушен поток; Protok zraka; Proud vzduchu; Онуов; Πορ αφο; Légáramlás; Aershreabhaidh; Gaisa plúsmá; Oro srautas; Fluss tal-arja; Debit de aer; Pnietok vzduchu; Pretok zraka.	Min.68.5l/min @ 5 PSI	Min.37.3l/min @ 5 PSI	Min.37.3l/min @ 5 PSI	Min.50.3l/min @ 5 PSI	Min.26.6l/min @ 1 PSI	Min.26.6l/min @ 1 PSI	Min.24.2l/min @ 1 PSI	Min.22.3l/min @ 1 PSI	Min.22.3l/min @ 1 PSI	Min.22.3l/min @ 1 PSI	Min.29l/min @ 1 PSI	Min.26l/min @ 1 PSI	Min.54.4l/min @ 1 PSI
Filter Efficiency; Efficacité du filtre; Filterleistung; Efficacia del filtro; Eficiencia del filtro; Filterefficiënte; Eficiência do filtro; Filteringsseffektivitet; Filteringsseffektivitet; Filter Effektivitet; Suodattimen tehokkuus; Výdajnosť filtra; Efectivitatea na filtrare; Učinkovitost filtra; Účinnost filtra; Filtri efektiivsus; Amfónómótio gátporo; Szűrőhatékonyság; Eifeachtacht Scagaire; Filtra eifeaktívlté; Filtr eifektivumas; Efficénza tal-filtra; Eficentja filtrului; Účinnost filtra; Učinkovitost filtra.	Min.99.99% @ 15L/min	Min.99.99% @ 30L/min	Min.99.99% @ 30L/min	Min.99.99% @ 30L/min	Min.99.99% @ 30L/min	Min.99.99% @ 30L/min	Min.99.99% @ 30L/min	Min.99.99% @ 30L/min	Min.99.99% @ 30L/min	Min.99.99% @ 30L/min	Min.99.99% @ 30L/min	Min.99.99% @ 30L/min	Min.99.995% @ 30L/min
Bacterial Efficiency; Rétenion bactérienne; Bakterienretention; Ritenzione di batteri; Retención de bacterias; Bacterieretentie; Retenção bacteriana; Bakterieretention; Bakterieretention; Bakterieretensjon; Bakteenien suodatus; Zatrzymywanie bakterii; Эффектност на бактериите; Učinkovitost bakterija; Bakteriální účinnost; Bakterite eifektivsus; Вострокуяі отомокоуэмкiтiтiо; Bakteriális hatékonyaság; Eifeachtacht Bactéaragh; Bakteriý eifektivlté; Bakteriý eifektivumas; Efficénza Batterjali; Eficentja bacterianá; Bakteriálna účinnost; Učinkovitost bakterij.	Min. 99.999%	Min. 99.999%	Min. 99.999%	Min. 99.999%	Min. 99.999%	Min. 99.999%	Min. 99.999%	Min. 99.999%	Min. 99.999%	Min. 99.999%	Min. 99.999%	Min. 99.999%	Min. 99.999%
Viral Efficiency; Rétenion virale; Virenretention; Ritenzione di virus; Retención de virus; Virusrelenie; Retenção de virus; Virusrelenie; Virusretention; Virusretensjon; Virusten suodatus; Zatrzymywanie wirusów; Вироуна ефективност; Virusna učinkovitost; Virusová účinnost; Virusi eifektivsus; Iwí отомотлоуэмкiтiтiо; Virushatékonyaság; Eifeachtacht Viressach; Virusu eifektivlté; Virusinis efektyvumas; Efficénza Viral; Eficentja viralá; Virusová účinnost; Virusna učinkovitost.	Min. 99.999%	Min. 99.999%	Min. 99.999%	Min. 99.999%	Min. 99.999%	Min. 99.999%	Min. 99.999%	Min. 99.999%	Min. 99.999%	Min. 99.999%	Min. 99.999%	Min. 99.999%	Min. 99.999%
Material Data; Matériau; Material; Dati materiale; Material; Material; Material; Materiale; Materiale; Material; Material; Дани на материала; Podaci o materijalu; Materiálové údaje; Materiaalsed andmed; Δεδομένα υλικοί; Anyagi adatok; Sonrai Ábhar; Matériale dati; Esminalni duomenys; Dejlja Materjali; Date materiale; Materiálové údaje; Materialni podatki.													
Housing; Boltier; Gehäuse; Alloggiamento; Carcasa; Behuizing; Estrutura; Kabinet; Hus; Hus; Pidke; Obudowa; Kópny; Kucište; Kryt; Úmbris; Периáлпуа; Burkolat; Tithiocht; Májokji; Bústas; Djar; Carcasă; Kryt; Ohisje.	PP	SBC	SBC + PVC + PA66	SBC	SBC	SBC	SBC	SBC	SBC	SBC	SBC	SBC	SBC
Filter Media; Matière du filtre : Filter-Medium : Materiale del filtro : Filtro : Filter medium : Matériau du filtre : Filtermaterial : Filtermateriale : Filtermedium : Suodattimen väline : Medium filtrage; Филтратуа среда; Filtrski medij; Filtračni prostředi; Filter; Μέσο φίλτρησης; Szűrőközeg; Meán Scagaire; Filtrét multividi; Filtruoti laikmena; Iifiltra iMidja; Medi de filtrare; Filtračné médium; Filtrimi medij.	Hydrophobic PTFE 1.0µm	Hydrophobic PTFE 1.0µm	Hydrophobic PTFE 1.0µm	Hydrophobic PTFE 1.0µm	Hydrophobic PTFE 1.0µm	Hydrophobic PTFE 1.0µm	Hydrophobic PTFE 1.0µm	Hydrophobic PTFE 1.0µm	Hydrophobic PTFE 1.0µm	Hydrophobic PTFE 1.0µm	Hydrophobic PTFE 1.0µm	Hydrophobic PTFE 1.0µm	Hydrophobic PTFE 1.0µm